AILORIA

DEPLAQUED

BEDIENUNGSANLEITUNG
INSTRUCTION MANUAL
MODE D'EMPLOI
ISTRUZIONI PER L'USO
MANUAL DE USUARIO
GEBRUIKSAANWIJZING
INSTRUKCJA OBSŁUGI

Version 1.0 - 10/2021

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
© VEHNS GROUP GmbH

PRODUCT FEATURES

The AILORIA DEPLAQUED sonic scaler is ideal for cleaning teeth at home. It removes tartar, plaque and discolouration, even from hard-to-reach areas.

- Intelligent high-frequency vibration
- Ergonomic design for comfortable handling
- 5 intensity levels
- Cordless battery operation
- Charging time approx. 2 hours
- Battery life approx. 6-8 hours when fully charged
- Charging of the lithium-ion battery via USB
- Rated voltage: 3.7V
- Rated current: 200mA
- Battery capacity: 300mAh
- Charging voltage: 5V

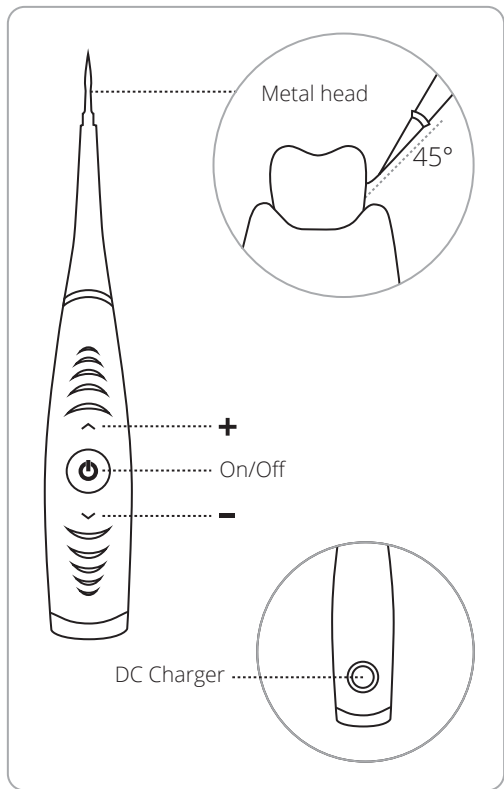SAFETY INSTRUCTIONS

- Always keep the USB cable and USB adapter away from water and other liquids.
- Remove the USB charging cable before cleaning the ultrasonic scaler.
- Do not use the product if it is damaged. The sonic scaler can only be repaired by the manufacturer.
- Children and persons with reduced physical or mental capabilities or lack of experience and knowledge should only use the product under adult supervision and instruction. Always keep DEPLAQUED out of the reach of children.
- DEPLAQUED is only suitable for cleaning teeth and should not be used for any other purpose.
- If excessive bleeding of the gums or pain occurs and persists while using the sonic scaler, please consult your dentist.
- The product should not affect pacemakers in any way. However, as a maximum precaution, we recommend that you consult your doctor before use.
- If you have undergone jaw or gum surgery in the last two months or if you have any health problems, consult your doctor before using DEPLAQUED.
- Clean the DEPLAQUED metal attachment with

water after each use and wipe the housing with a damp cloth. Sterilize the tooth cleaner with alcohol if necessary.

- The sonic scaler is protected against water splashing according to IPX6 protection standard. However, do not keep the product under water for a long period of time.
- The product contains a rechargeable battery. Keep the sonic scaler away from fire and high temperature to avoid danger.

COMMISSIONING

Charge the sonic scaler before the first use.

Charging time is approximately 2 hours to charge the sonic scaler before first use or fully discharged battery.

CHARGING PROCEDURE

The connector for the charging cable is located at the bottom of the back of the unit. Open the silicone cap and connect the charging cable to the sonic scaler. Connect the other end to a power source (wall socket) using the USB charging adapter provided. Alternatively, you can also connect the USB cable to a USB interface (e.g. PC or notebook).

The sonic scaler cannot be used during the charging process. A red light shines in the transparent base

for the duration of the charging process. When the charging process is completed, this is indicated by a green light.

APPLICATION

Press the ON/OFF button to switch on the sonic scaler. Select the desired intensity using the arrow keys.

The following 5 intensity levels are available:

Intensity level 1: 8321 vibrations per minute

Intensity level 2: 9787 vibrations per minute

Intensity level 3: 11195 vibrations per minute

Intensity level 4: 12423 vibrations per minute

Intensity level 5: 13320 vibrations per minute

Start cleaning the teeth by gently placing the metal attachment on the areas with plaque and discolouration. Now slowly move the metal tip over the areas to be cleaned. This process can be repeated several times until the desired result is achieved. The sonic scaler can also be used on the edges of the teeth and in the spaces between the teeth. Move the sonic scaler carefully to prevent possible injuries to the gums.

Clean the metal head with water after each use and wipe the housing with a damp cloth. Sterilize the sonic scaler with alcohol if necessary.

DISPOSAL

The product must not be disposed of with normal household waste. Follow local regulations for separate disposal of electrical and electronic products.

Proper disposal will help prevent negative effects on the environment.

MANUFACTURER

VEHNS GROUP GmbH
Theatinerstrasse 40-42
80333 Munich, Germany

service@vehnsgroup.com
www.vehnsgroup.com
www.ailoria.de

MADE IN CHINA

Subject to alterations and errors

The manufacturer accepts no liability for damage caused by incorrect, improper or incompatible use of the product.

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VEHNS GROUP GmbH, Theatinerstrasse 40-42